

## Photomultiplier Base with Preamplifier

- For use with 10-stage PMTs that fit standard 14-pin sockets
- Built-in low-noise preamplifier
- Both preamplifier output and anode output
- Test input for system testing
- Protection circuit for internal transistors
- Internal gain adjustment



The ORTEC Model 276 Photomultiplier Tube Base and Preamplifier incorporates an integral low-noise preamplifier, a PMT base with voltage divider network, and a focus control for optimum performance in scintillation detector applications. The unit provides two outputs: the preamplifier output for energy analysis and the anode output for either timing or auxiliary energy analysis. The preamplifier is dc-coupled to simplify pole-zero cancellation in the main amplifier. A Test input accepts the output of a pulse generator to calibrate and test the preamplifier and the following system. The Model 276 has a diode protection network to prevent damage to the internal

transistors due to sudden application or removal of high voltage to the unit. A simple internal modification in the unit allows the gain to be adjusted for any value desired by the user. The Model 276 is powered from any ORTEC main amplifier or preamplifier power supply.

The Model 276 is directly compatible with 10-stage PMTs that fit standard 14-pin sockets including those listed in Table 1.

The Model 276 is also compatible with other 10-stage tubes not listed in Table 1 (see Fig. 1). Compatibility may be determined by comparison with those listed.

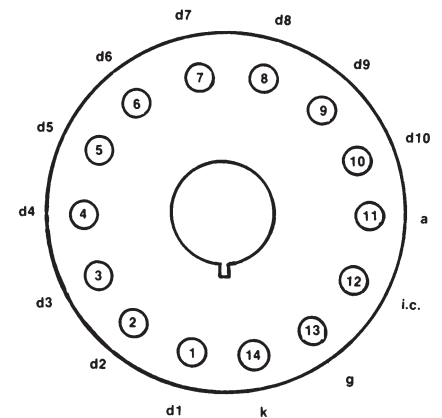


Fig. 1. JEDEC B14-38 PMT Pin Base, with Pin Assignments:

- d1 – d10 dynodes 1 to 10
- a anode
- i.c. internal connection
- g grid
- k cathode

Table 1. Compatible Photomultiplier Tubes.

ADIT		Burle (formerly RCA)		Hamamatsu		Philips
B51B01	B76C01	4900	S83020F	PM55	R1513	XP2202
L51B01	B89B01	5819	S83021E	R208	R1612	XP2203B
V51B01	B89C01	6342A	S83022F	R550	R1791	XP2412B
B51D01	B89D01	6655A	S83025F	R594	R1836	
B51C01	B133D01	S83006E		R877	R1847-07	
B76B01	B133C01	S83013F		R878	R1848-07	
V76B01	V133B01	S83019F		R1507	7696	
				R1512		

## Specifications

### PERFORMANCE

#### PREAMPLIFIER

**Integral Nonlinearity**  $<\pm 0.02\%$ , 0 to +10 V.

**Temperature Instability**  $<\pm 0.005\%/^{\circ}\text{C}$ , 0 to  $50^{\circ}\text{C}$ .

**Output Rise Time**  $<100$  ns for test input or fast scintillator.

**Output Fall Time** Time constant of 50  $\mu\text{s}$ .

**Output Noise**  $<50$   $\mu\text{V}$  rms with ORTEC Model 572 Amplifier and time constant of 1  $\mu\text{s}$ .

**Conversion Gain** Nominally 5  $\mu\text{V}/\text{eV}$  with 2- by 2-in. NaI(Tl) crystal and PMT gain of  $10^6$ ; the typical output for a 511-keV x ray with a 10-stage PMT gain of  $10^5$  will be  $\sim 250$  mV.

**Saturation Level** +10 V into an open circuit; +5 V into 93- $\Omega$  load.

**VOLTAGE DIVIDER** Resistor-divider connected to 10-stage PMT base. Total resistance 1.49 M $\Omega$ , resulting in bleeder current of 0.6 mA with typical high voltage of 1 kV. The distribution is linear to all stages with the focus adjustment on the grid.

### CONTROL

**FOCUS** Single-turn locking potentiometer on panel for external adjustment of PMT grid potential.

### INPUTS

**POS HV** SHV connector, AMP 51494-2, for distribution of positive high voltage to PMT base; +2000 V maximum.

**TEST** BNC connector, accepts pulses from an ORTEC pulse generator for testing and calibration.

**SIGNAL** Preamplifier input is connected internally to dynode 10.

**POWER** Captive 3-m (10-ft) power cable terminated in Amphenol 17-20090 connector accepts preamplifier operating power; compatible with all ORTEC main amplifiers and the Model 4002P Portable Power Supply.

**PM SOCKET** TRW 3B14. Fits JEDEC B14-38 PMT pin base (see Fig. 1).

### OUTPUTS

**PREAMP** BNC connector, furnishes preamplifier positive output pulse to an ORTEC main shaping amplifier for linear energy analysis;  $Z_o = 93$   $\Omega$ , dc-coupled.

**ANODE** BNC connector, furnishes negative anode output pulse for use either in timing or auxiliary energy analysis;  $Z_o = 1$  k $\Omega$ , ac-coupled.

### ELECTRICAL AND MECHANICAL

**POWER REQUIRED** For preamplifier, +24 V, 16 mA; -24 V, 16 mA; for PMT base, +2000 V maximum (use rated voltage for the tube that is installed).

### WEIGHT

**Net** 0.65 kg (1.5 lb).

**Shipping** 1.3 kg (3.0 lb).

**DIMENSIONS** 5.6 cm (2.2 in.) diam x 10.2 cm (4 in.) long plus 3 m (10 ft) captive power cable.

## Accessories

A C-36-12 Cable is required for connection to a high-voltage supply. A C-24-12 Cable is recommended for connecting the preamplifier output to the spectroscopy amplifier. For timing signals from the anode, use a C-25-12 Cable.

## Ordering Information

To order the Model 276 Photomultiplier Base with Preamplifier or related accessories, use the following model numbers and descriptions:

Model	Description
276	Photomultiplier Base with Preamplifier
C-24-12	RG-62A/U 93- $\Omega$ Cable with two BNC male plugs, 12-ft length
C-25-12	RG-58A/U 50- $\Omega$ Cable with two BNC male plugs, 12-ft length
C-36-12	RG-59A/U 75- $\Omega$ Cable with two SHV female plugs, 12-ft length

Specifications subject to change  
072108

**ORTEC**<sup>®</sup>

[www.ortec-online.com](http://www.ortec-online.com)

Tel. (865) 482-4411 • Fax (865) 483-0396 • [ortec.info@ametek.com](mailto:ortec.info@ametek.com)  
801 South Illinois Ave., Oak Ridge, TN 37831-0895 U.S.A.  
For International Office Locations, Visit Our Website

**AMETEK**<sup>®</sup>  
ADVANCED MEASUREMENT  
TECHNOLOGY